



---

# AI Derma Checker: Skin Disease Detection Platform

User Mannual

---

## Contents

<b>1</b>	<b>Introduction</b>	<b>2</b>
<b>2</b>	<b>System Requirements</b>	<b>2</b>
<b>3</b>	<b>Accessing Derma AI</b>	<b>2</b>
<b>4</b>	<b>How to Use Derma AI: Step-by-Step Guide</b>	<b>2</b>
4.1	Step 1: Landing on the Website . . . . .	2
4.2	Step 2: Uploading an Image . . . . .	3
4.3	Step 3: Payment Process . . . . .	4
4.4	Step 5: AI Analysis and Viewing Results . . . . .	5
<b>5</b>	<b>Frequently Asked Questions</b>	<b>6</b>
<b>6</b>	<b>Troubleshooting</b>	<b>6</b>
<b>7</b>	<b>Support and Contact</b>	<b>7</b>
<b>8</b>	<b>Legal Disclaimer</b>	<b>7</b>

# 1 Introduction

Welcome to **Derma AI** – a smart dermatology assistant powered by Deep Learning. Upload a clear image of a skin condition, make a secure payment, and receive a fast AI-generated prediction for possible skin diseases.

**Disclaimer:** This tool is for informational purposes only. Always consult a certified dermatologist for a medical diagnosis.

# 2 System Requirements

- Active internet connection
- Modern web browser (Chrome, Firefox, Safari, Edge)
- Valid image file (.jpg, .jpeg, .png)
- Payment method (Credit/Debit Card)

# 3 Accessing Derma AI

1. Open a browser and go to: <https://www.dermaai.com>
2. You'll arrive at the homepage with a brief overview.

# 4 How to Use Derma AI: Step-by-Step Guide

This guide covers the complete user flow from uploading an image to downloading your diagnostic report in five clear steps.

## 4.1 Step 1: Landing on the Website

- Visit the Derma AI application homepage.
- Click the **Upload** button at the top right to begin.

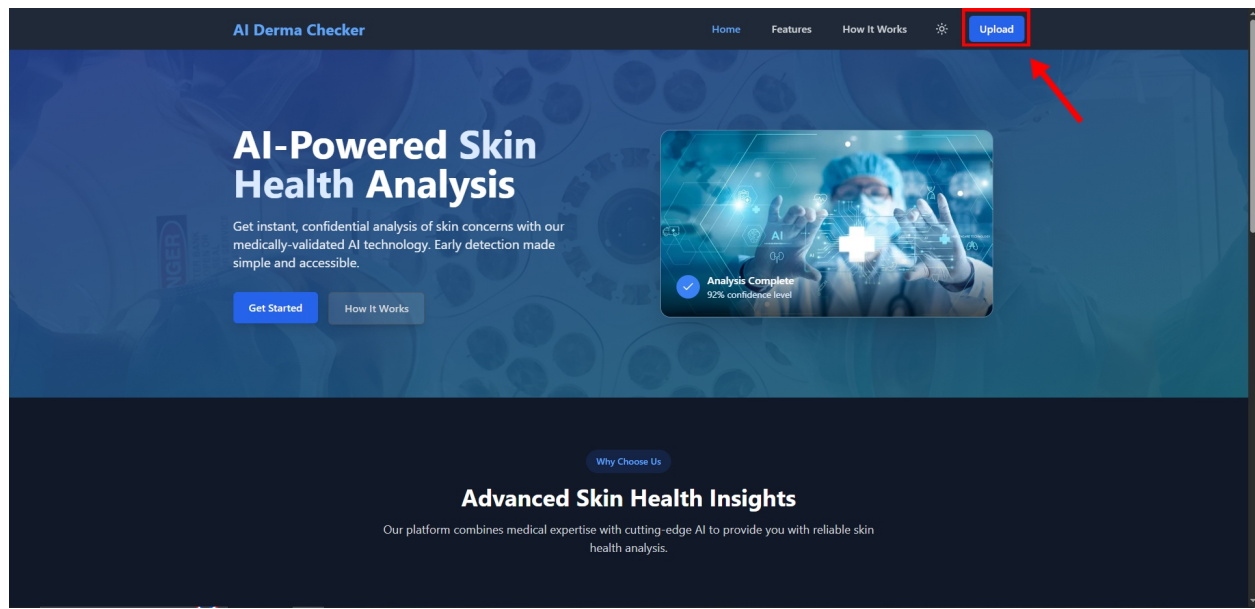


Figure 1: Homepage with Upload Button

## 4.2 Step 2: Uploading an Image

1. Choose a clear image of the affected skin area.
2. Review the image preview:
  - If the image is fine, click **Proceed to Payment**.
  - If not, click **Cancel** and re-upload a better image.

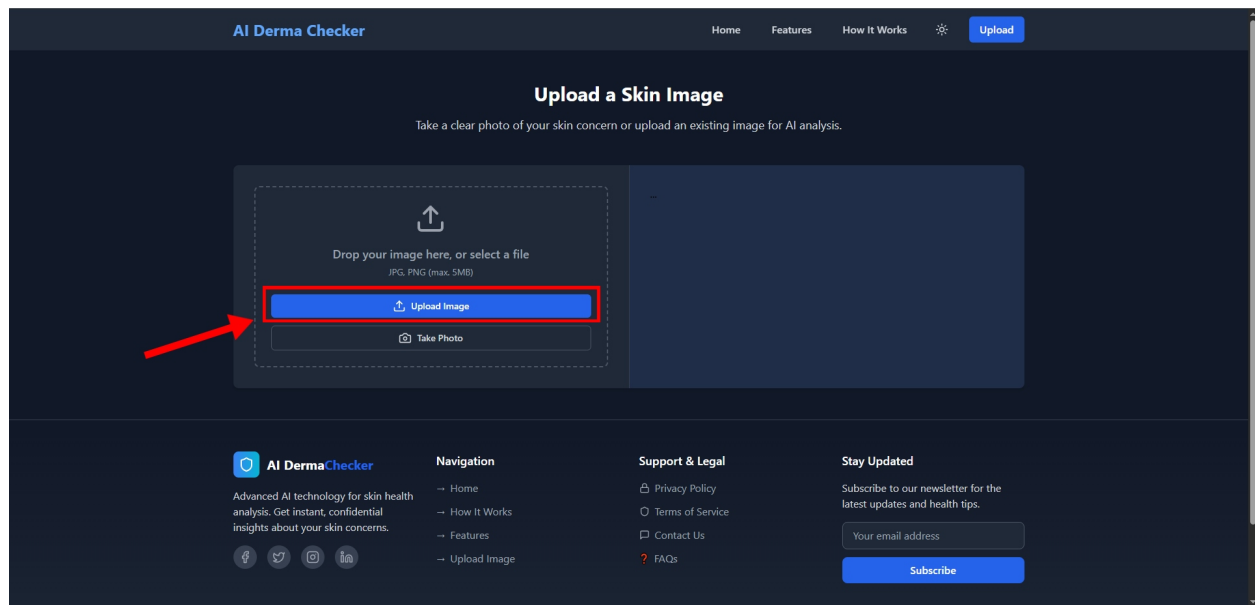


Figure 2: Image Upload and Preview

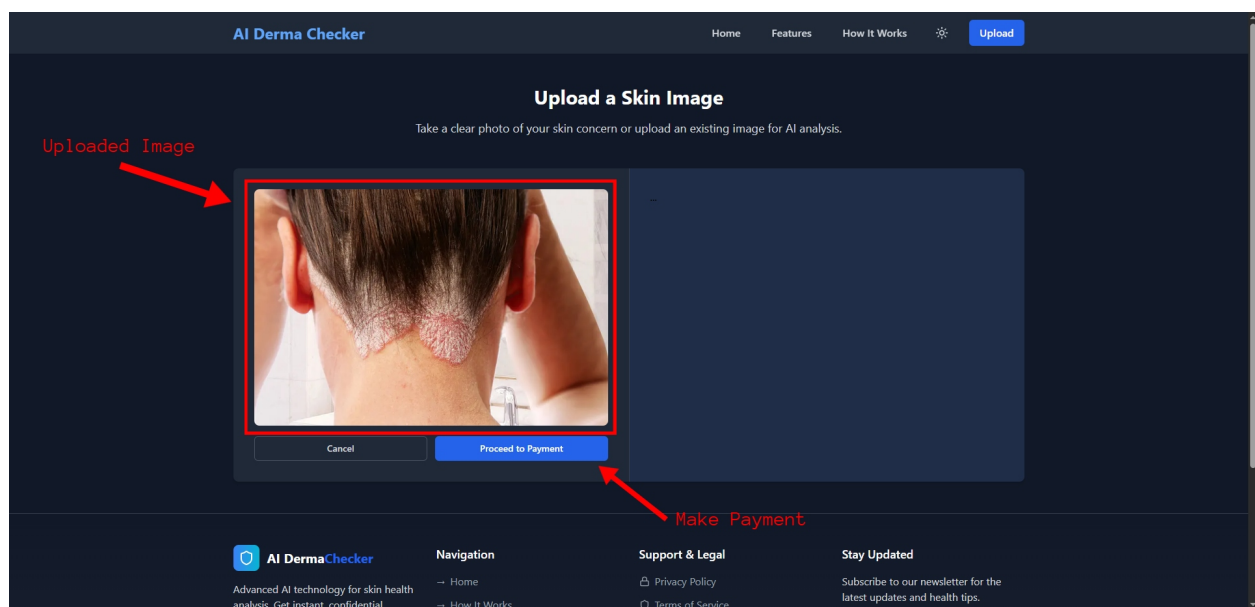


Figure 3: Image Upload and Preview

*Tip: Ensure the image is in focus, well-lit, and shows the area clearly.*

### 4.3 Step 3: Payment Process

1. You'll be redirected to a secure Stripe payment page.
2. Enter valid card details and confirm the payment.

### 3. Scenarios:

- **Success:** You'll be redirected to the results page.
- **Invalid Details:** Stripe displays an error, allowing retry.
- **Insufficient Funds:** Stripe denies the payment.
- **Payment Gateway Error:** You may retry later.
- **Browser Closed:** The session is terminated, and no result is generated.

The screenshot displays the Stripe Secure Payment Page. On the left, there's a sidebar with a 'New business sandbox' header and a 'Choose a currency' dropdown set to 'PKR 1,467.84'. Below this, it shows '1 USD = 293.5680 PKR' and 'Skin Disease Detection' with 'PKR 1,467.84'. The main area on the right features a 'link' logo with '4444', an 'Email' field, a 'Card Information' section with a card number, cardholder name 'Abdullah Butt', and country 'Pakistan'. Below this is a checkbox for 'Securely save my information for 1-click checkout' and a 'Pay' button. Red arrows point to the 'Email' field and the 'Pay' button, labeled 'User Details' and 'Make payment' respectively.

Figure 4: Stripe Secure Payment Page

*Note: All transactions are secured via Stripe. Payment is required to proceed.*

## 4.4 Step 5: AI Analysis and Viewing Results

- Once payment is successful, AI processing is triggered.
- The system takes around 5–15 seconds to analyze the image.
- then, You can view the results on-screen.
- Pdf can be download which consists of Predicted skin condition, Confidence percentage and Medical etc.
- Click **Save in PDF** to save a copy of your diagnosis report.

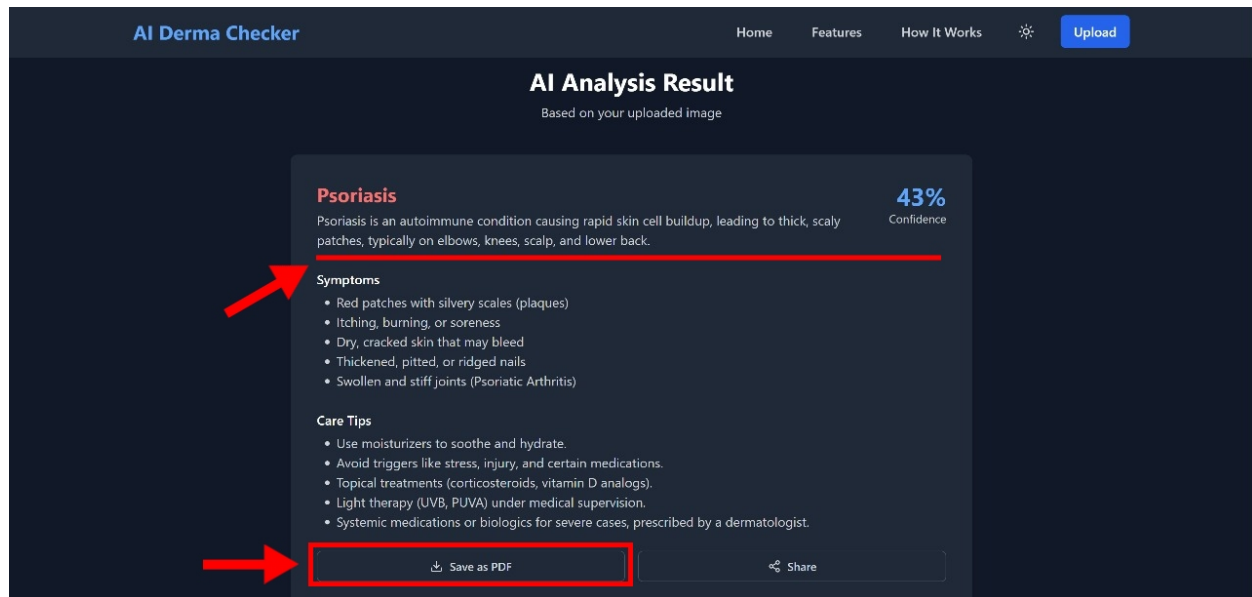


Figure 5: Diagnosis Results with Download Option

*Note: The PDF report contains your result and treatment recommendation for future reference.*

## 5 Frequently Asked Questions

### Q1: How accurate is Derma AI?

Derma AI is trained on a large and diverse dermatological dataset. While it provides high-confidence predictions, it is not a substitute for a professional medical opinion.

### Q2: Are my images safe and private?

Yes. Your images are encrypted and transmitted securely. We do not store or share any personal data.

### Q3: Can I get a refund?

Refunds are only offered in case of technical errors or no diagnostic output.

## 6 Troubleshooting

- **Image not uploading:** Check internet connection and supported file format.
- **Payment failure:** Verify card details or try a different method.

- **No result displayed:** Refresh the page and try again. If unresolved, contact support.

## 7 Support and Contact

For any issues, questions, or feedback:

- **Email:** support@dermaai.com
- **Phone:** +1-234-567-890
- **Live Chat:** Available on our website

## 8 Legal Disclaimer

Derma AI is an informational tool and does not replace professional medical advice, diagnosis, or treatment. Always seek the advice of a qualified healthcare provider with any questions regarding a medical condition.